

# Technical Data Sheet

Doc No–Ab/TDS/BC-147

Rev No: 02 Date 02/2022-23

## TECHNICAL DATA SHEET

**Abdos Labtech Pvt. Ltd.**

**Aspirator Bottle with Stopcock, PP**

**Valid for Catalog No. P11515, P11516, P11517**

### **1. Description/Specification**

#### **1.1 Description:**

Aspirator Bottle with Stopcock, PP

1.1.1 Aspirator bottles are manufactured using high purity virgin USP Class VI Medical grade polypropylene (PP).

1.1.2 These bottles come with secure sealing cap for better security and long-term storage.

1.1.3 These bottles are used for storing and shipping reagents, chemicals, culture media, and standards.

1.1.4 These are very strong and durable for extended use over time with a strong carrying handle.

1.1.5 Comes with high quality and durable stopcock or tap for easy withdrawal of liquids, Chemicals, reagents, culture media and standards.

1.1.6 Can be used at a temperature below 0°C.

1.1.7 Superior chemical resistance

#### **1.2 Raw material:**

**1.2.1 Body:** PP (Poly Propylene)

**1.2.2 Cap:** HDPE (High Density poly Ethylene)

#### **1.3 Material:**

1.3.1 Poly Propylene (PP), Hydrophobic, confirm US FDA 21CFR & High Density poly Ethylene (HDPE), Hydrophobic, confirm US FDA 21CFR

1.3.2 Complies with the requirement of European Pharmacopoeia, Plastic Containers and closures for Pharmaceutical use.

1.3.3 The product has passed the USP testing including class VI tests ,Materials has successful passed the biological tests according to ISO 10993 external communication devices for indirect blood contact for Prolonged period.

1.3.4 Free from Heavy metals, RoHS Free according to meets relevant requirements of EU regulations 1935/2004/CE ,2002/95/EC as amended, UE 10/2011 as amended. EU CONEG & it's following amendments. It is safe use in contact with food stuff, Pharmaceuticals & drinking water

1.3.5 Product is manufactured without the use of Raw Materials of animal origin. Therefore BSE/TSE issue does not concern this product to the best of our Knowledge. No materials containing latex or natural rubber are used in the Manufacturing, handling & packaging process for this product.

#### **1.4 Quality Control:**



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

1.4.1 Product are manufactured, inspected and found to be in compliance with product & quality specification requirements as documented in ISO 9001:2015 & 13485:2016 Quality Management System. Testing of various characteristics in accordance with the valid specification, inspection reviewed & signed for product release.

1.4.2 Can be used at a temperature below 0°C.

## 1.5 DIMENSIONAL SPECIFICATIONS:

Code	Capacity (Ltr)	Mouth Diameter (mm)	Height (mm)	No. per Case
P11515	5	51.4	374	1
P11516	10	51.4	445	1
P11517	20	51.4	525	1

## 1.6 GENERAL FEATURES:

CE	Yes
Autoclavability	Yes
Translucence	Yes
Resistant	To be effected by radiations
Sterilization	NA
Temperature	From 0°C to 121°C
Leak Proof	Yes

1.7 Chemical resistance: See website ([www.abdoslifesciences.com](http://www.abdoslifesciences.com))

1.7.1 Manufactured in Aseptic Environment

1.8 Other Information (not applicable if neglected):

1.8.1 Fair resistance; exposure (30 days) to the chemical does not cause any damage.

1.8.2 Sufficient resistance; exposure causes damage of poor importance, which sometimes is only temporary.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

**1.8.3 Poor resistance;** exposure to chemical is not allowed, or causes immediate permanent damage.

Results of testing with reference materials are to be considered as indication: in case of specific use it is recommended to carry out preliminary testing.

Information relating to the resistances of plastic materials to temperatures, sterilization and washing treatments, have been made as per reports of the raw materials' producers, following the data published in the literature and the experience gained in the use of products.

## **WARNING**

**DO NOT** store strong oxidizing agents.

**DO NOT** place any plastic labware in a flame.

**DO NOT** mix any chemicals that may result in a thermic reaction, which can cause product failure

Consult local fire codes prior to storage of flammable liquids in carboys/bottles/aspirators.

Consult OSHA guidelines prior to handling or lifting filled carboys/bottles/aspirators.

**If applicable:** Using the spigot as a handle can ruin the threads and causing the spigot not to seal properly.

## **Cleaning:**

Bottles/Aspirators/Carboys need to be cleaned by hand using a non-abrasive, neutral pH, mild detergent that does not contain a sheeting agent.

## **Cap and Spigot Cleaning:** (if applicable)

1. Remove O-ring from cap.
2. Disassemble the spigot.
3. Spray and wipe the spigot, cap and O-ring with 70% Isopropyl Alcohol.
4. Rinse all parts with distilled water.
5. Replace O-ring in cap.
6. Autoclave cap and spigot at 121°C at 15 psi for 20 minutes.

## **Instruction for Use**

Product can be used for food contact also.

Storage temp. Max 30°C

Before using the product check crackness, pouch punctured etc.

Check all the information printed on pouches or boxes.

Always wear gloves while handling specimen.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakheswari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA


**Registered office**  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903  
**Corporate Office**  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7

# Technical Data Sheet

## UKRP

Obelis UK LTD  
Sandford Gate,  
East Point Business Park,  
Oxford OX4 6LB  
United Kingdom  
Tel: +44,1491.378012  
info@obelis.co.uk

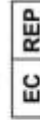

Obelis s.a.  
Bd General Wahis 53  
1030 Brussels Belgium  
Tel: +(32) 2 732-59-54  
Fax: +(32) 2 732-60-03  
mail@obelis.net



## Abdos Labtech Pvt. Ltd.

Plot No: 1 Lakshwari Bhagwanpur Roorkee 247667 (UK) India  
Email : labtech@abdosindia.com, www.abdoslifesciences.com  
Ph: +(91) 11 41081897



Quality Assurance Incharge  
Abdos Labtech Pvt. Ltd.



Abdos Labtech Private Limited  
(A Division of Abdos Group)  
(An ISO 9001,13485,14001,18385,GMP & CE certified company)  
Shiv Ganga Industrial Estate  
Plot No-1, Lakshwari, Bhagwanpur,  
Roorkee-247667, (UK) INDIA

Registered office  
13A Government  
Place East, 2nd Floor Kolkata  
West Bengal – 700069  
CIN: U85100WB2008PTC123903

Corporate Office  
319, 3rd Floor, DLF Tower A, DDA Dist. Centre  
Jasola, New Delhi - 110025

Ph. No: +91- (0) 1332235277  
Email:labtech@abdosindia.com  
Web: [www.abdoslifesciences.com](http://www.abdoslifesciences.com)  
GST: 05AAGCA8605A1Z7